


CHINESE SMALL VOLUME (CSV) ACCESSORIES FOR VISION® G2 TESTERS

Teledyne Hanson Chinese Small Volume (CSV) components adapt directly to Vision® G2 dissolution testers for applications requiring less than 500 mL of dissolution media and to conform with the specifications of the Chinese Pharmacopoeia (CP). The CP-compliant 250 mL CSV vessel has a greater capacity than the traditional 150 mL small-volume vessel. Two ports in the vessel cover allow manual sampling or auto-sampling through a fixed vessel-mount probe and an optional return-to-vessel port. When configured with the Hanson SuperMag™, autosampling can also be performed through magazine-mounted auto-probes.

CHINESE SMALL VOLUME WITH MANUAL OR AUTOMATED SAMPLING

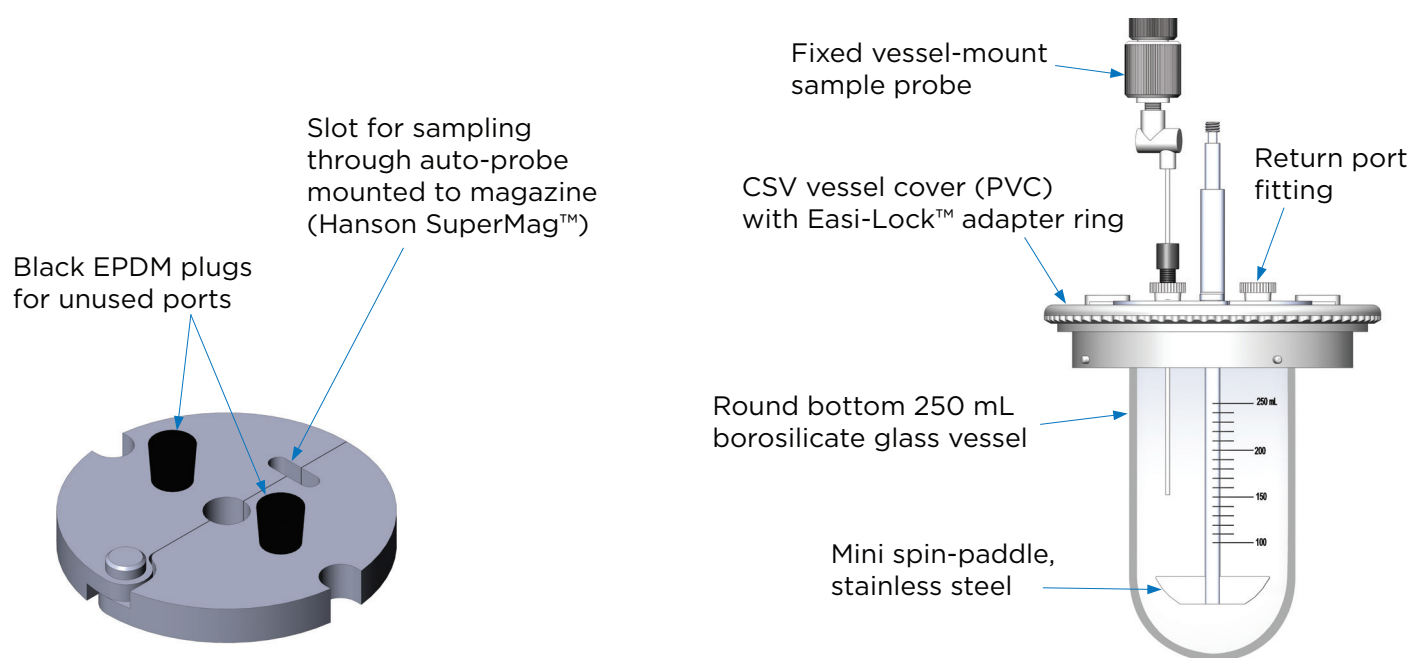

ORDERING INFO

For further information please contact your Teledyne Hanson Research representative or email us at hansonsales@teledyne.com.

CSV VESSELS AND COVERS

- 74-105-280** CSV Vessel, Round Bottom, 250 mL, Clear, Vision G2
- 74-105-281** CSV Vessel, Round Bottom, 250 mL, Amber, Vision G2
- 74-104-173** CSV Vessel Cover, CS, Vision G2 (Compatible with SuperMag™)
- 74-104-175** CSV Adapter Ring Assembly, Vision G2 (required for 250 mL Chinese Small Volume)

CSV PADDLES AND PROBES

- 74-105-290** CSV Mini Spin-Paddle, 316 SS, Vision G2
- 74-104-201** Manual Sampling Cannula, Adjustable, 1/8 in., PEEK
- 74-104-211** Fixed Vessel-Mount Sample Probe, Adjustable, 1/16 in., 316 SS (requires Tubing Harness)
- 74-104-205** Return Port Fitting
- 74-105-307** Auto-Probe (1), 1/8 in., PEEK (requires SuperMag™ and Tubing Harness)