


MEDIA REPLACEMENT FOR VISION G2[®] CLASSIC 6[™]

MEDIA REPLACE BEAKER SYSTEM

The Vision[®] G2 Classic 6[™] is an award-winning 6-position dissolution tester. With the addition of the new Media Replace Beaker system, the Classic 6 can now simultaneously warm an additional 1 liter of media while performing a six-sample dissolution test. The warmed media can then be replaced either manually into the dissolution vessels or automatically using the Vision[®] AutoPlus[™] autosampler.

The **Media Replace Beaker System** consists of a jacketed beaker that warms media to the intended test media temperature. The two ports in the beaker's jacket connect it to both the Classic 6 waterbath and Vision Heater, which allows warmed water to continuously flow through the entire system, creating a stabilized temperature environment for media.

The Media Replace Beaker system is an optional accessory available for your new or existing Classic 6 dissolution tester.

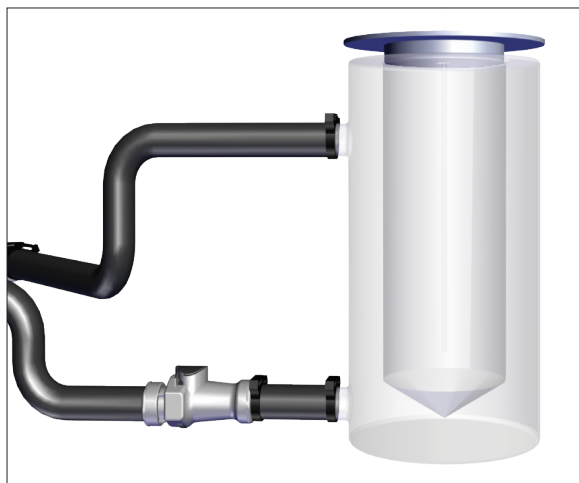

FEATURES

- Adds a seventh position for warmed media
- Jacketed exterior keeps media warm for immediate use during your test
- Comes complete with a cover and 6 probes to connect to your Vision[®] AutoPlus[™]
- Includes all the tubing and quick-disconnect fittings to connect to your Classic 6[™]
- Field retrofit option available
- Holds 1 L of solution

ORDERING INFO

For more information, please contact your local Teledyne Hanson representative or email us at:

hansonsales@teledyne.com

74-207-015 Media Replace Beaker System, Classic 6



9810 Variel Avenue
Chatsworth, CA 91311 USA

Phone +1 818.882.7266
www.teledynehanson.com

© Teledyne Hanson, a division of Teledyne Instruments, Inc. | 99-310-054 Rev. 03 January 2020